


MORPHEUS®
SMART PICC

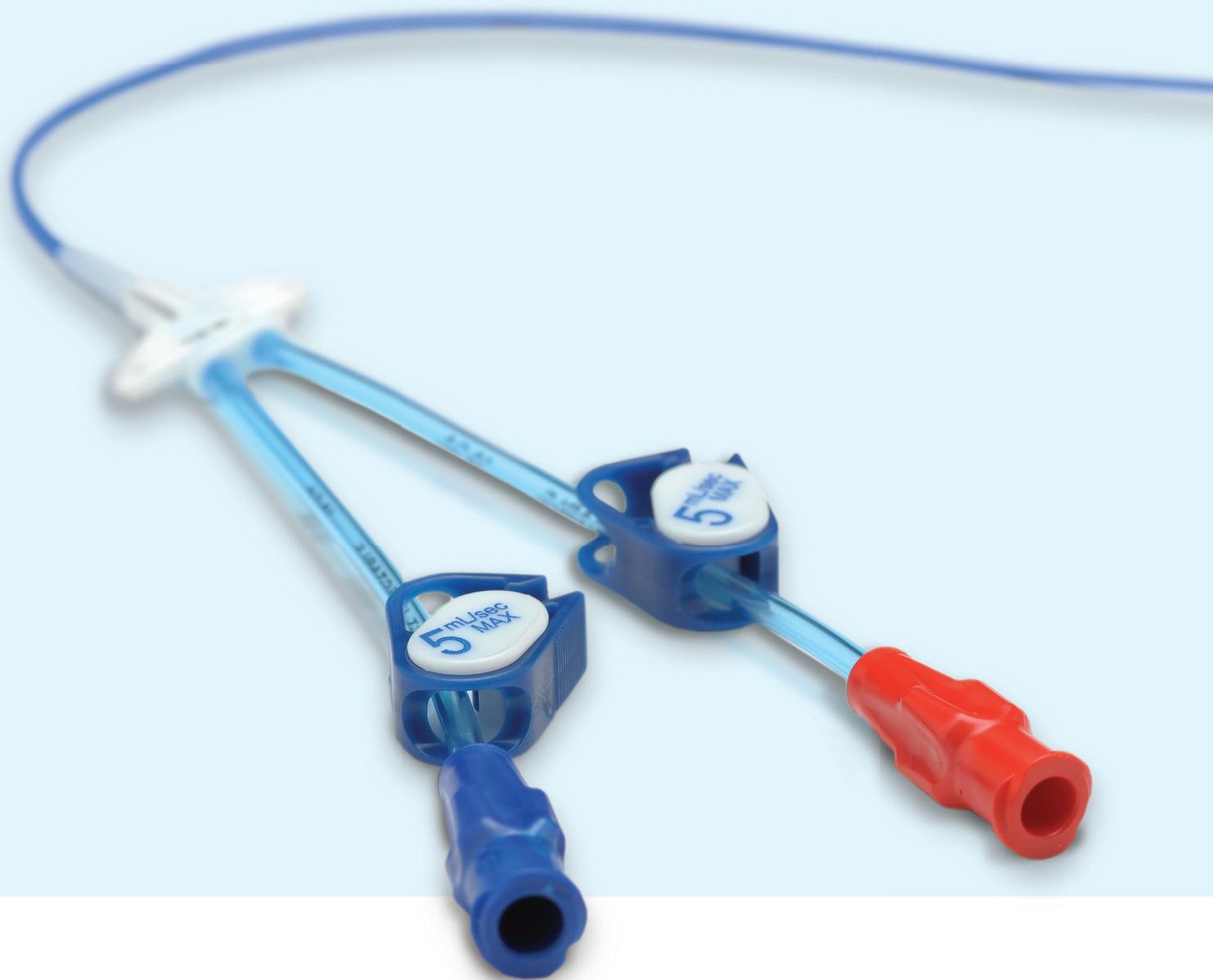

ANGIODYNAMICS®

Engineered for Life™


The Morpheus® SMART PICC provides clinicians with a comprehensive line of CT compatible PICC kits that offer True French size catheters with the shortest taper in the industry; this results in less trauma to the vessel, increased blood flow to the vessel and ultimately decreases the patient's risk of thrombosis, infection and device occlusion.

Best in Class Flow Rate
CT Rated and labeled for CT injection.



All-Carbothane® catheter
The chemically resistant all-Carbothane catheter offers proven drug compatibility.

6F Catheter Comparison

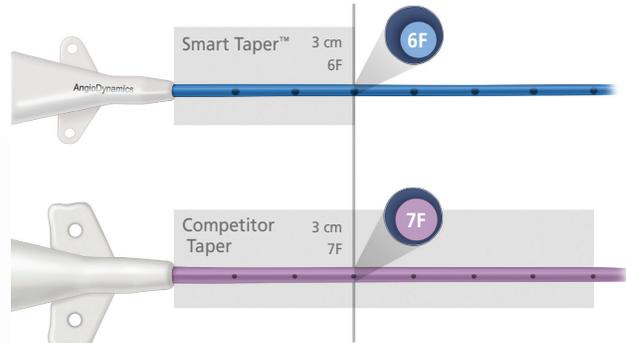

Smart Taper™ Technology

The Smart Taper™ reverse taper technology is designed to reduce the risk of thrombosis and vessel trauma.


Highly Visible Material


Superior product visibility, complete with radiopaque markings with centimeter spacing provide increased visibility, easier catheter placement and tip confirmation.

True French Sizing

True French size catheter means less catheter in the patient and allows for an atraumatic PICC insertion.

5F Dual Catheter Profile

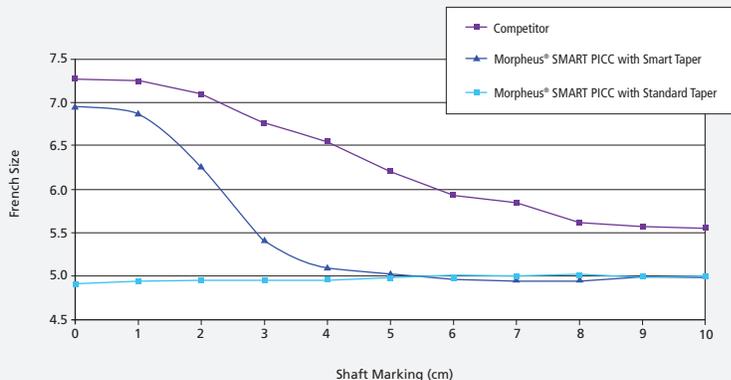

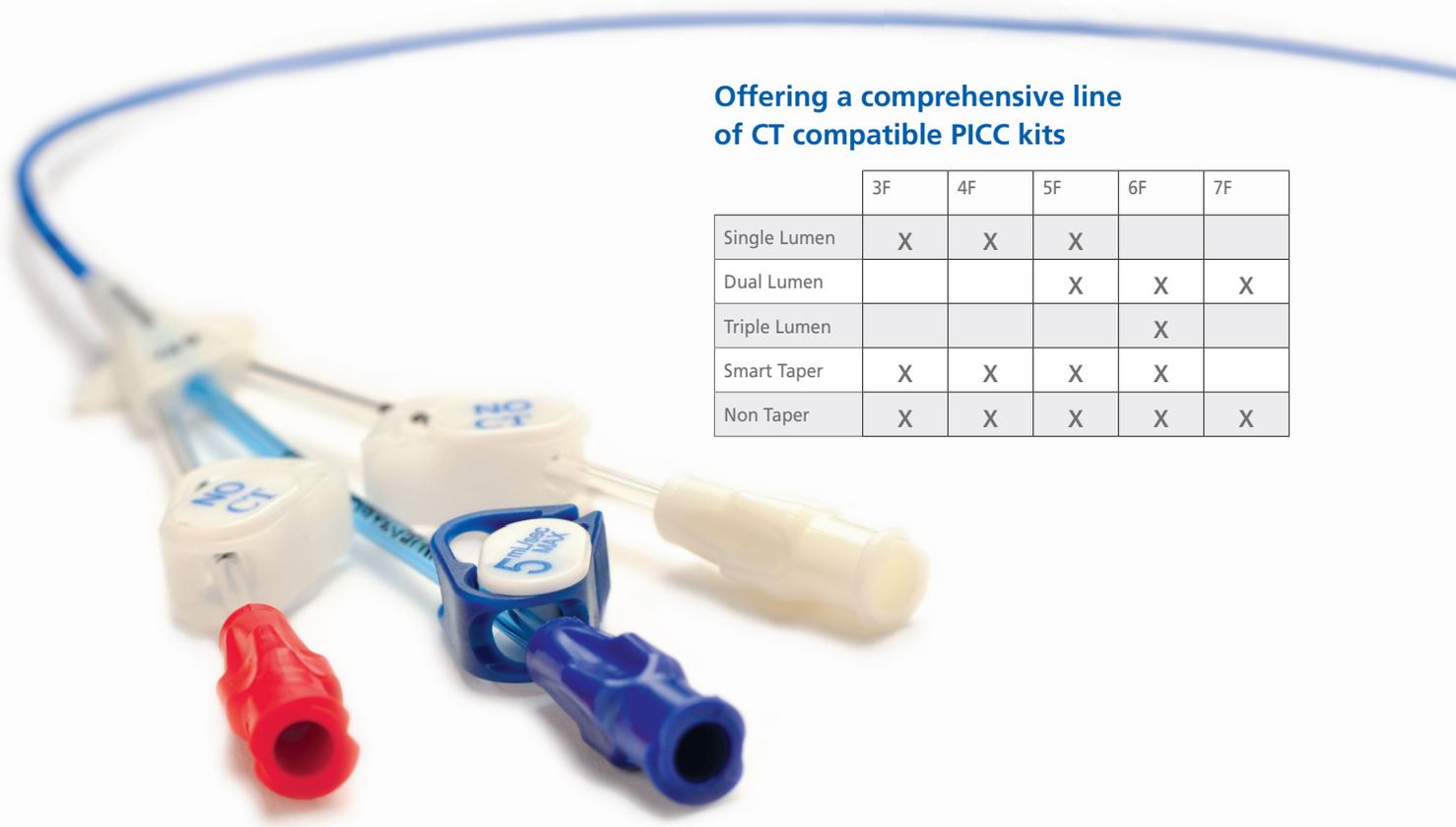

Offering a comprehensive line of CT compatible PICC kits

	3F	4F	5F	6F	7F
Single Lumen	X	X	X		
Dual Lumen			X	X	X
Triple Lumen				X	
Smart Taper	X	X	X	X	
Non Taper	X	X	X	X	X

IMPORTANT RISK INFORMATION

INDICATIONS FOR USE:

The AngioDynamics Morpheus® SMART PICC CT is indicated for short or long term peripheral access to the central venous system for intravenous therapy, power injections of contrast media, and allows for central venous pressure monitoring. For blood sampling, infusion, or therapy use a 4 French or larger catheter. For central venous pressure monitoring, it is recommended that a catheter lumen of 20 gauge or larger be used.

CONTRAINDICATIONS:

This device is contraindicated whenever the patient

has or is suspected of having: infection, bacteremia, or septicemia; allergic reactions to the device's materials; insufficient body size to accommodate the size of the implanted device; past irradiation; previous episodes of venous thrombosis or vascular surgical procedures at the prospective placement site; local tissue factors will prevent proper device stabilization and/or access.

WARNINGS AND PRECAUTIONS:

Please see package insert for complete list of warnings and precautions.

POTENTIAL COMPLICATIONS:

Please see package insert for complete list of potential complications.

Indications, contraindications, warnings and instructions for use can be found in the instructions for use supplied with each device. Observe all instructions prior to use. Failure to do so may result in patient complications.

CAUTION: Federal (USA) law restricts these devices to sale by or on the order of a physician.

ANGIODYNAMICS®

USA > 603 Queensbury Avenue, Queensbury, NY 12804 > tel: 800-772-6446 or 518-798-1215 > fax: 518-798-1360
 International > Building 2000, Beach Drive, IQ Cambridge, Waterbeach, Cambridge, CB25 9TE > United Kingdom
 tel: +44 1223 729 300 > fax: +44 1223 729 349

www.angiodynamics.com

AngioDynamics and Morpheus are registered trademarks of AngioDynamics, Inc. Smart Taper is a trademark of AngioDynamics, Inc. All other trademarks or registered trademarks are the property of their respective owners. U.S. patents pending.



Proudly distributed by:

MEDICAL SPECIALTIES AUSTRALIA PTY LIMITED

54 GIBBES STREET, CHATSWOOD NSW 2067
 PHONE 02 9417 7955 FAX 02 9417 5779
 E-Mail: custserv@msa.com.au
www.msa.com.au

ADC1068-01/2013